


MonoLead ECG Lead-Wire Set



**More time with your
patients, less time
with wires**

“I need an ECG lead-wire set that’s
easy to apply, won’t tangle, and is
comfortable for my patients.”

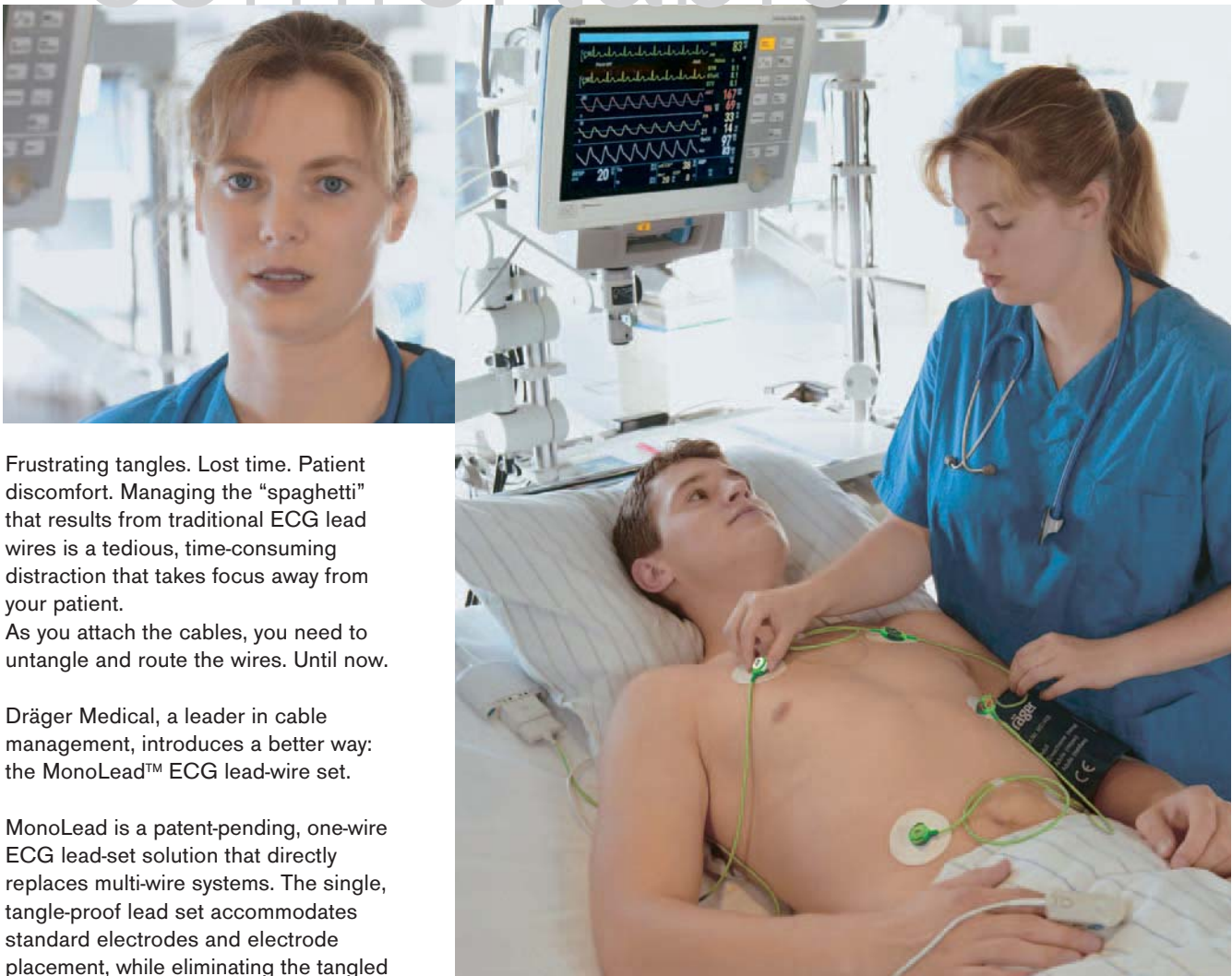

Frustrating tangles. Lost time. Patient discomfort. Managing the “spaghetti” that results from traditional ECG lead wires is a tedious, time-consuming distraction that takes focus away from your patient.

As you attach the cables, you need to untangle and route the wires. Until now.

Dräger Medical, a leader in cable management, introduces a better way: the MonoLead™ ECG lead-wire set.

MonoLead is a patent-pending, one-wire ECG lead-set solution that directly replaces multi-wire systems. The single, tangle-proof lead set accommodates standard electrodes and electrode placement, while eliminating the tangled clutter that frequently occurs with traditional multi-lead sets. Ideal for both adult and pediatric patients, MonoLead is compatible with all Dräger Medical Infinity® patient monitors.



More comfortable for your patients, easier for you



Because MonoLead is a single cable rather than a complicated clutter of wires, it's less intimidating for your patient and their loved ones. Additionally it is more physically comfortable. That's because MonoLead is designed with Z-Snap™ zero-insertion force electrode snaps, which connect easily and hold securely, with no pressure or discomfort for your patient.

MonoLead is fast and easy to apply, and uses the standard lead placement that's familiar to nurses everywhere. It's flexible, durable and easy to clean. And the bright green color makes the cable easy to identify – even under a sheet.

So spend less time with wires and more time with your patients. Discover MonoLead, the one-wire ECG cable solution that's better for you, better for your patients.

Only MonoLead offers these benefits:

- Tangle proof – replaces 3-, 5- and 6 electrode lead sets with a single wire
- Fast and easy to apply – has clearly labeled, color-coded snaps
- Same placement and results – supports the standard (Mason-Likar) electrode placement on the torso that is familiar to clinicians worldwide
- More comfortable – incorporates a patented Z-Snap™ zero-insertion force electrode snaps
- Suitable for both pediatric and adult patients
- Easy to clean – wipes clean with standard disinfectants

www.monolead.com

Ordering Information

Pods:		3368391; 5950196	5191221	MS14688	MS16162	5589663
		MultiMed 5, 2,5m	MultiMed 6	MULTIMED 6T	MULTIMED 6T OR	MultiMed 12
		Multimed 5, 1,5m				
MonoLeads T-connector						
MS14556	MonoLead 3 ECG CBL T-connector IEC2	x (a)	x (a)	x	x	
MS14560	MonoLead 5 (limb) ECG CBL T-connector IEC2	x (a)	x (a)	x	x	
MS14682	MonoLead 6 ECG CBL T-connector IEC2		x (a)	x	x	
Adapters:						
MS14679	Adapter MonoLead to MultiMed 5	x				
MS14680	Adapter MonoLead to MultiMed 6		x			
HF-protection for OR						
5947226	ESU block for Multimed 5	x				
7486140	ESU block for Multimed 6		x			
MonoLead Dualpin						
MS16233	MonoLead 3 ECG CBL D/PIN IEC2	x	x			
MS16229	MonoLead 5 (limb) ECG CBL D/PIN IEC2	x	x			x
MS16230	MonoLead 5 (chest) ECG CBL D/PIN IEC2					x

x (a) Adapter required

USA:

Draeger Medical, Inc.

3135 Quarry Road
Telford, PA 18969

USA

Tel: (215) 721-5400

Toll-free: (800) 437-2437

Fax: (215) 723-5935

E-mail: wwwinfo@draegermed.com

Manufacturer:

Draeger Medical, Inc.
Telford, PA 18969, USA

The quality management system at
Draeger Medical, Inc. is certified according
to ISO 13485, ISO 9001 and Annex II of
Directive 93/42/EEC (Medical devices).

www.draegermedical.com