

## Smart Solutions Work Place – Endotracheal Suctioning

Planned process-support is a major contribution to time and cost savings within today's medical environment. Allow us to help you, through our global experience, to ensure that the work processes around your patients are the best available.



### SMART SUCTIONING

- Ergonomic positioning – everything within easy reach
- Standardized setup – harmonized workflows
- Clear separation – improved hygiene
- Secure device mounting - help prevent tripping and spillage

### SUCTIONING SET-UP

Mounted approximately 950 mm from the floor, a shelf with side rails supports a sterile working surface. The vacuum regulator is mounted on the shelf-rail, on the patient side, the display clearly visible. In front of the regulator, also on the rail, a catheter holder places the catheter tip within easy reach. Wire baskets for catheter- and sterile glove storage can be positioned within arms length, on the shelf rail.

The secretion jars are mounted on a standard rail attached to a swiveling arm beneath the shelf. During suctioning or when the jars have to be emptied, the arm can be extended forward, positioning the secretion jars at the patient bedside. When not in use they can easily be swung back under the work surface, saving valuable work space and providing hygienic separation from the patient.

The swiveling arm and rail can also be used for positioning other accessories that have the same functional requirements, such as a litter bag or sharps bin.

## ORDER INFORMATION

### Parts required

Parts required	Ordernumber
Support arm for front rail	G 96 526
or	
support arm for 38 mm tube	G 96 527
Rail holder, short	G 13 822
Individual standard rail, 600 mm	G 13 829
Small equipment pole, 500 mm	G 15 676
Clamp ring set	G 19 015
Bronchial aspiration for system VarioVac® B/OP rail	MP 00 505
Vacuum tube corresponding with national standards	86 01 697
Bronchial aspiration system VarioAir® B/OP rail	MP 00535
Medical compressed air hose corresponding with national standards	86 01 697
Secretion jar set	2M 85 056
Ceiling mount	2M 85 119
Aspiration hose holder	2M 85 239
Rail holder	M 26 842
<b>Additional options</b>	
Connecting-hose holder	2M 85 446
Wire basket for disposable gloves, 200 mm	M 26 146
Wire basket for catheter storage, 480 mm	M 25 121
Holder for litter bags incl. 100 bags	M 24695




MT-172-2006

Work Place –  
Endotracheal Suctioning



D-28534-2009

VarioAir® B/OP rail

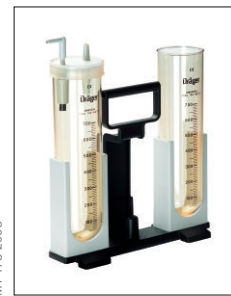

MT-173-2003

Secretion jar set



2-568-95

Wire basket



MT-0007-2004

VarioVac® vacuum regulator



MT-3713-2008

Holder for litter bag

### HEADQUARTERS

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany

[www.draeger.com](http://www.draeger.com)

### REGION EUROPE CENTRAL AND EUROPE NORTH

Dräger Medical GmbH  
Moislinger Allee 53–55  
23558 Lübeck, Germany  
Tel +49 451 882 0  
Fax +49 451 882 2080  
[info@draeger.com](mailto:info@draeger.com)

### REGION EUROPE SOUTH

Dräger Médical S.A.S.  
Parc de Haute  
Technologie d'Antony 2  
25, rue Georges Besse  
92182 Antony Cedex, France  
Tel +33 1 46 11 56 00  
Fax +33 1 40 96 97 20  
[dlimfr-contact@draeger.com](mailto:dlimfr-contact@draeger.com)

### REGION MIDDLE EAST, AFRICA, CENTRAL AND SOUTH AMERICA

Dräger Medical GmbH  
Branch Office Dubai  
Dubai Healthcare City, P.O. Box 505108  
Dubai, United Arab Emirates  
Tel + 971 436 24 762  
Fax + 971 436 24 761  
[contactuae@draeger.com](mailto:contactuae@draeger.com)

### REGION ASIA / PACIFIC

Dräger Medical  
South East Asia Pte Ltd  
25 International Business Park  
#04-27/29 German Centre  
Singapore 609916, Singapore  
Tel +65 6572 4388  
Fax +65 6572 4399  
[asia.pacific@draeger.com](mailto:asia.pacific@draeger.com)

### REGION NORTH AMERICA

Draeger Medical, Inc.  
3135 Quarry Road  
Telford, PA 18969-1042, USA  
Tel +1 215 721 5400  
Toll-free +1 800 437 2437  
Fax +1 215 723 5935  
[info.usa@draeger.com](mailto:info.usa@draeger.com)

### Manufacturer:

Dräger Medical GmbH  
23542 Lübeck, Germany  
The quality management system at  
Dräger Medical GmbH is certified  
according to ISO 13485, ISO 9001  
and Annex II.3 of Directive 93/42/EEC  
(Medical devices).