


## Dräger DCD 5000 Oral fluid collection system

Simple, fast and secure detection of intoxicants and drugs:  
The Dräger DCD 5000 does not only enable an easy collection of saliva (oral fluid), you can also store the samples and send them to the laboratory for analysis – without the sample being falsified along the way.

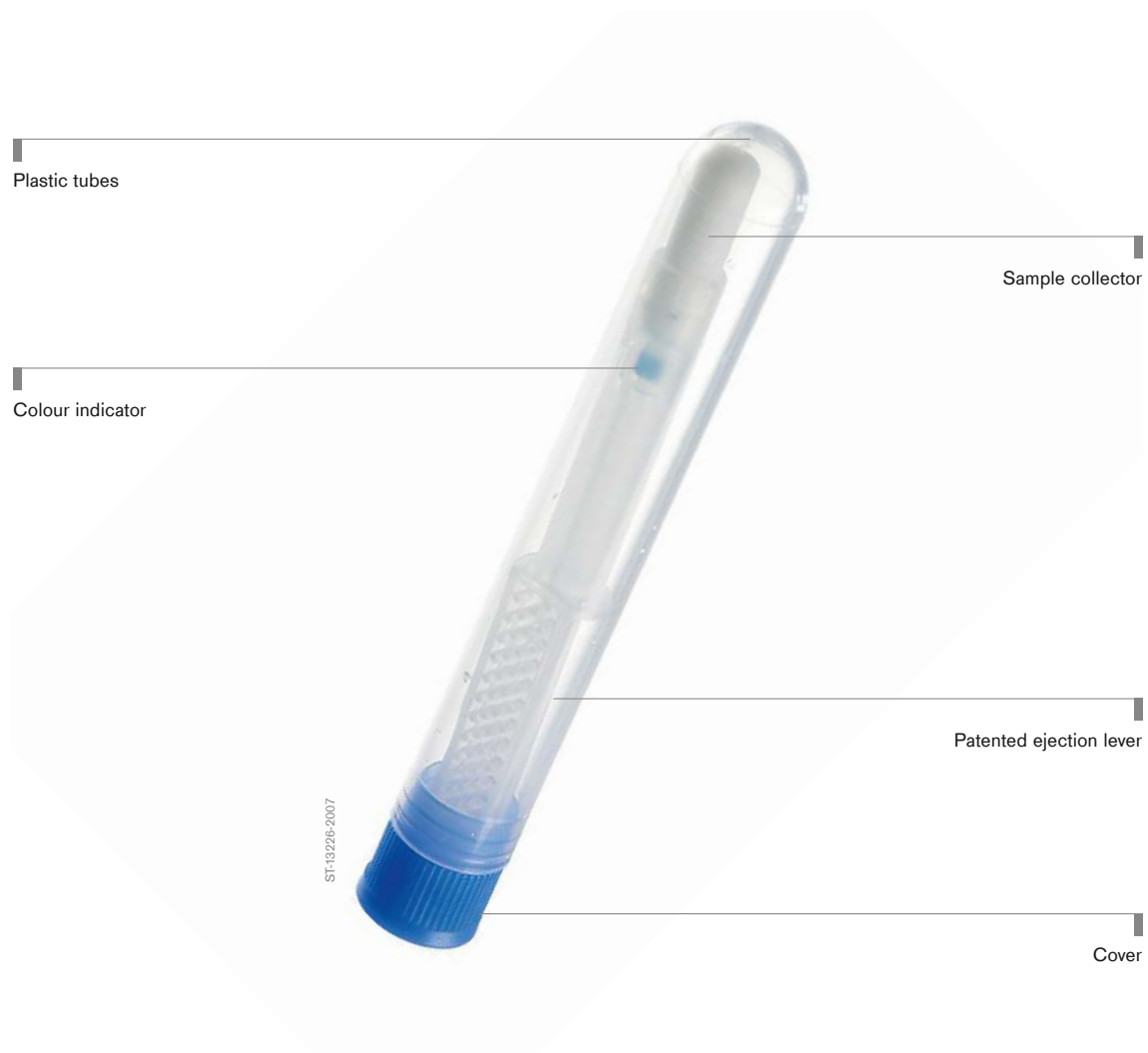

## Benefits

---

### **Ultra filtrate of blood**

Saliva is an established ultra filtrate in clinical chemical diagnostics. Saliva samples are approved as sample material for the detection and diagnostic analysis of substance abuse.

---

### **Check without interruptions**

The continuous check of sampling is facilitated with the method based on saliva, compared to urine sampling. Thereby, manipulation by the test person is eliminated.

---

### **Easy, clean and safe**

The results of a pretest are confirmed reliably via a laboratory analysis. Storage, transport and recovery of the sample must not impair this result. The DCD 5000 sampling kit allows for this security.

---

### **Double sample**

Combining two DCD 5000 samplers, it is very easy to collect two separate samples simultaneously. Afterwards, the samples can be separated and serviced separately, e.g. for transport and storage.

## Details



D-15392-2017

Remove the sampler from the tube and give it to the test subject.



D-15404-2017

To collect the saliva sample, swab from one inner cheek to the other.



D-15353-2017

When the colour indicator turns blue, the sampling can be ended.



ST-13225-2007

The sample is secured inside the tube for further use.



ST-13230-2007

Patented ejection mechanism to release the mouthpiece from the handle.



ST-13228-2007

Twin sampler to collect two saliva samples at once.

## Accessories



D-15335-2017

### Dräger SCSK

In the SCSK confirmation set, the DCD 5000 is enhanced with a stabilising component (isopropanol) to ensure high-accuracy detection of drugs in the sample, even after lengthy periods between sampling and laboratory testing.

## Accessories



D-15278-2017

### Dräger SCK

The SCK confirmation set contains two samplers and stabilising components (isopropanol) as well as documentation forms and shipping materials, ensuring a chain of custody through clear documentation and secure shipment of samples to the lab.

## Related Products



D-42702-2015

### Dräger DrugCheck 3000

Use the Dräger DrugCheck 3000 to find out within minutes if a person recently consumed certain drugs. The compact and quick saliva (oral fluid) based drug test yields reliable results affordably and easily. The non-invasive drug test kit does not require electricity and can therefore be used anywhere.



D-54720-2012

### Dräger DrugTest 5000

No pipetting, no drips, no timing: With the Dräger DrugTest 5000, a drug test is carried out simply and quickly. The collected saliva sample can be analyzed immediately – for accurate results on the spot.

## Technical Data

Dimensions of the individual sampler (H x W x D)	120 mm x 15 mm x 15 mm
Working range	Storing an unused Dräger DCD 5000: 4 to 30 °C
Sample volume	approx. 0.4 ml of saliva (oral fluid)
Licenses/Standards	For non-medical purposes or forensic applications (USA only)

## Ordering Information

	<b>Packaging unit</b>	<b>Part no.</b>
Dräger DCD 5000 Sampler Kit	20 saliva samplers, packaged individually in plastic tubes with seal, 1 x instructions for use	83 19 910
SCSK – Saliva Confirmation Single Kit	Set consists of 10 sealed plastic bags, each containing 1 x DCD 5000 saliva sampler, 1 x tube of isopropanol, 1 x instructions for use included per set	83 26 300
SCK – Saliva Confirmation Kit	Set consists of 2 x DCD 5000 saliva samplers, 2 x tubes of isopropanol with packaging material: 1 x sealable transport bag for diagnostic material, 2 x safety seals, 1 x carbon copy form for laboratory use and instructions for use	83 23 880

## Notes

Not all products, features, or services are for sale in all countries.  
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to [www.draeger.com/trademarks](http://www.draeger.com/trademarks) to find the current status.

**CORPORATE HEADQUARTERS**  
Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany  
[www.draeger.com](http://www.draeger.com)

**REGION EUROPE**  
Dräger Safety AG & Co. KGaA  
Revalstraße 1  
23560 Lübeck, Germany  
Tel +49 451 882 0  
Fax +49 451 882 2080  
[info@draeger.com](mailto:info@draeger.com)

**REGION MIDDLE EAST, AFRICA**  
Dräger Safety AG & Co. KGaA  
Branch Office  
P.O. Box 505108  
Dubai, United Arab Emirates  
Tel +971 4 4294 600  
Fax +971 4 4294 699  
[contactuae@draeger.com](mailto:contactuae@draeger.com)

**REGION ASIA PACIFIC**  
Dräger Singapore Pte. Ltd.  
61 Science Park Road  
The Galen #04-01  
Singapore 117525  
Tel: +65 6872 9288  
Fax: +65 6259 0398  
[asia.pacific@draeger.com](mailto:asia.pacific@draeger.com)

**REGION CENTRAL  
AND SOUTH AMERICA**  
Dräger Indústria e Comércio Ltda.  
Al. Pucurui - 51 - Tamboré  
06406-100 - Barueri - SP  
Tel. +55 (11) 4689-4900  
[relacionamento@draeger.com](mailto:relacionamento@draeger.com)

Locate your Regional Sales  
Representative at:  
[www.draeger.com/contact](http://www.draeger.com/contact)

