


Linea Living Elegant wall supply unit

With its timeless design and versatile lighting options, Linea Living is the all-purpose solution for lighting and media supply in the patient's room.

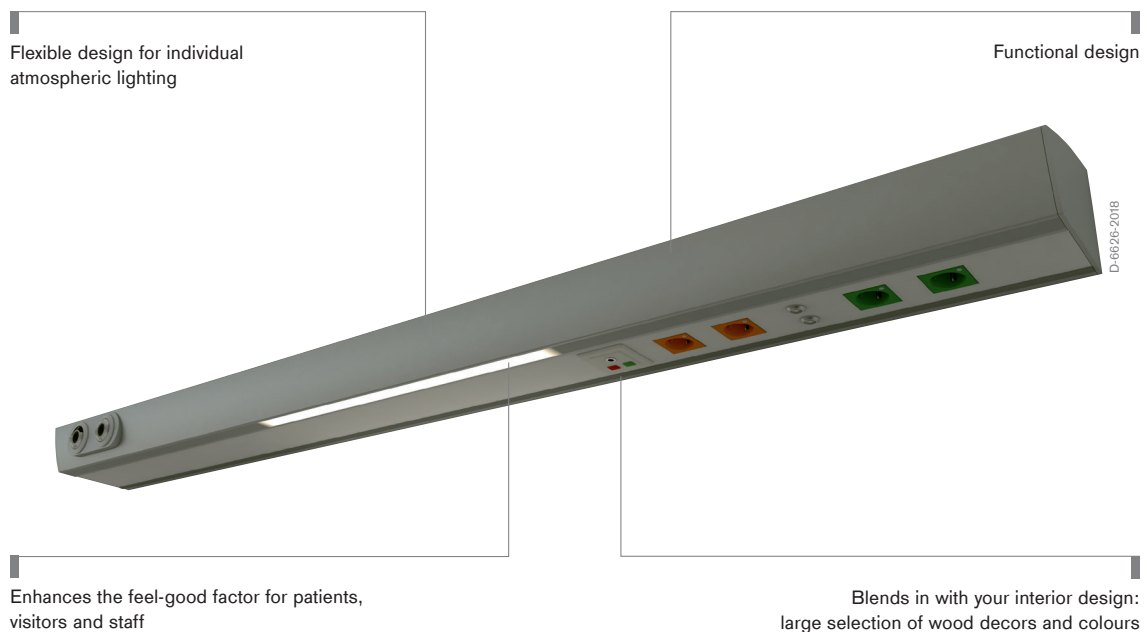

Product benefits

Solutions for you and your patients

Design your daily work processes more efficiently while also boosting the feel-good factor for your patients, visitors and staff. That is the goal of Linea Living and of our entire range of high-performance wall supply units. The products in this range can be flexibly adapted to your existing infrastructure and are highly compatible with one another. Their high degree of versatility gives you the freedom to create designs from purely functional to design-optimised.

Flexible lighting concepts

Linea Living provides glare-free light with high illuminance and offers excellent colour rendering and light distribution. A variety of light options (indirect room lights, reading lights or night lights) are available in various versions (colour temperature and light intensity) for the comfort and safety of both patients and staff.

For a wide range of applications

Linea Living can be used in many different ways. Such as an independent and comprehensive medical supply unit in care and rehabilitation wards or as an additional lighting unit, for example in combination with a Linea Vertica in an intensive care unit. The modular design of Linea Living ensures it is always the first choice.

Sleek, functional design

Linea Living consists of robust extruded aluminium profiles with rounded edges and perfectly integrated outlets for maximum safety, ease of cleaning and hygiene. At the same time it offers a range of outlet options, providing good access to gas, electrical and data terminals. Optionally, an equipment rail can be added for attaching medical equipment. Thanks to the free choice of different lengths, Linea Living is suitable for single bed as well as multiple-bed configurations.

Wide range of design options

Colours can enhance the well-being of patients and thereby promote the healing process whilst also creating a positive working environment for clinical staff. The housings and equipment panels are available in 17 standard RAL colours as well as various modern wood finishes. The slim, unobtrusive profile of Linea Living enables it to blend in with existing interior designs and colour schemes.

Easy installation and maintenance

The smart installation design makes installation and connection easier than ever. That makes Linea Living particularly well-suited for upgrading existing installations. Its innovative construction and modern LED technology make it very low-maintenance.

Details



D-7623-2018


Linea Living in combination with Linea Vertica



D-7624-2018

Feel-good atmospheres with countless design options

Colour and design options

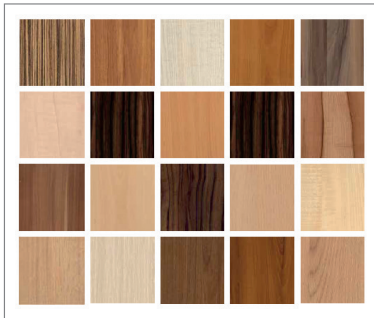

D-4121-2016


Selection of available wood decors



D-4120-2016

Selection of available RAL colours

Profile cross section

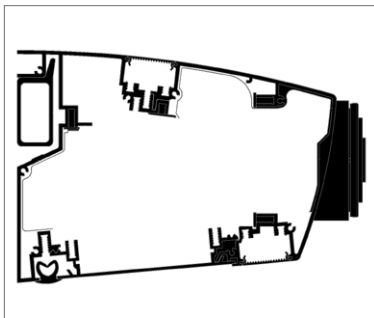

D-7627-2018

Profile cross section

Examples of use



Linea Living for illumination

When used as a dedicated illumination unit, Linea Living is an effective complement to other supply units such as Linea Vertica, Linea Ambience or even third-party systems. It can be fitted with a room light at the top, a reading light at the bottom plus a night light.



Linea Living for illumination with electrical outlets


In addition to the illumination options, the unit can be configured with both electrical and data outlets. The outlets can be freely positioned on the service plate on the underside of Linea Living in order to bring the services as close as possible to the patient. The lighting options available comprise of a room light, a reading light and a night light.




Linea Living as a full supply unit

As a full supply unit version, Linea Living, illumination, electrical and data outlets can be freely configured. In addition, terminal units can be positioned on the front of Linea Living. The underside is equipped with electrical and data outlets. As lighting solutions, a room light, a reading light and night light are available.

Light distribution

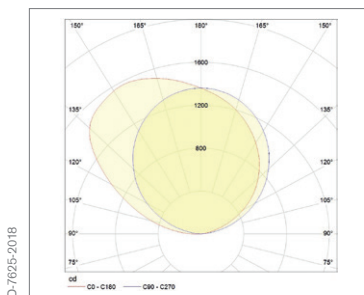


Light distribution curve of indirect lighting

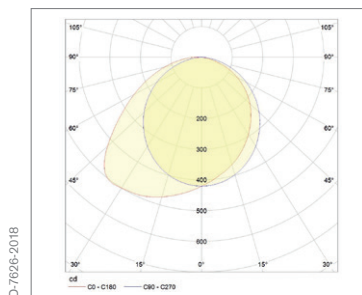

Light distribution curve of the reading light

Verwandte Produkte

D-11430-2016

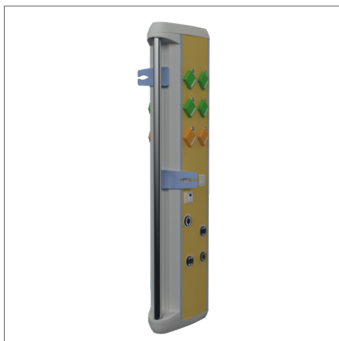

GeminaDUO

Whether for an emergency room, a recovery room or an intensive care unit, the Dräger GeminaDUO is just perfect for your workplace. The wall supply unit lets you choose from a wide range of configurations and options, helping to create an efficient and highly flexible workplace.

D-11431-2016



Linea

Linea N, Linea IM and Linea I offer the perfect combination of ergonomics and economy. These wall supply units can be individually configured to provide customised electrical, gas and IT supply solutions with illumination options for a variety of general and specialised applications.


D-18384-2017



Linea Vertica

Combining function and design – the Linea Vertica functional wall supply unit provides a high degree of flexibility for individual workstation design. The modern look and broad range of design and illumination options enable a patient-friendly atmosphere to be created – even with the most demanding clinical requirements.

D-18385-2017

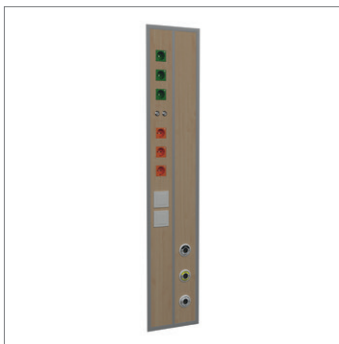

Linea Flush

Linea Flush is both modular and space-saving, so it can be used to meet the special requirements of the most diverse range of areas. The space-saving design enables it to be used for a wide range of applications beyond or in addition to conventional wall supply units.

Technical Data

Linea Living	
Height / depth	130 mm / 179 mm without terminal units), 195 mm (with terminal units)
Length	Your choice of any length (up to 6 m in a single piece)
Illumination	Indirect room light (switch or dimmer-controlled): warm white (3,000 K) with 5,400 lm, or warm white (3,000 K) with 7,600 lm, or neutral white (4,000 K) with 5,700 lm, or neutral white (4,000 K) with 8,000 lm Reading light (switch or dimmer-controlled): warm white (3,000 K) with 1,650 lm, or neutral white (4,000 K) with 1,750 lm Night light on the top of the Linea unit
Accessories	Dräger accessories for equipment rails (only for versions with optional equipment rail)
Rows of outlets	1 x below for electronic and communications installations / 1 x in front for terminal units
Components	Outlets in various colours, versions and national standards, potential equalisation pins, terminal units for medical gases, low-current equipment, such as nurse call, RJ45 data ports, telephone sockets, etc.
Wall mount	max. 65 kg/m, shear force per stud bolt (M10) 500 N.
Media connection	Type of installation and feed points as required – depending on the structural conditions
Colour options	17 standard RAL colours or various wood finishes
Load	The optional integrated equipment rail can support a maximum of 12.5 kg of additional load per rail.
Medical device	Complies with MDD (Medical Device Directive), medical electrical devices IEC 60601-1 and 60601-1-2, Medical Supply Units ISO 11197

Notes

Notes

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

MANUFACTURER:

V.T.S. Vision Technology Systems GmbH
Arnsberg, Deutschland
The quality management system of V.T.S Vision Technology Systems GmbH is certified according to ISO 13485, ISO 9001 and Annex II of Directive 93/42/EEC (medical devices).

REGION DACH

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA PACIFIC

Draeger Singapore Pte. Ltd.
25 International Business Park
#04-20/21 German Centre
Singapore 609916
Tel +65 6308 9400
Fax +65 6308 9401
asia.pacific@draeger.com

REGION CENTRAL AND SOUTH AMERICA

Draeger Panamá Comercial
S. de R.L.
59 East Street, Nuevo Paitilla,
House 31, San Francisco Town
Panama City, Panama
Tel +507 377-9100 / Fax -9130
servicioalcliente@draeger.com

Locate your Regional Sales Representative at:
www.draeger.com/contact

