

Dräger-Vapor® 2000 and D-Vapor® vaporizers for volatile anesthetic agents

With over 400,000 units sold around the world, Dräger vaporizers have set high standards for more than 50 years. Consistent, long-lasting quality made in Germany – by Dräger.



ECONOMY MEETS ERGONOMICS

Both the Vapor 2000 and the D-Vapor have been designed to accommodate 300 ml of liquid anesthetic. That is more than the entire volume of a standard bottle. You can even add 250 ml if the Vapor is not quite empty. This means you can easily refill a Vapor and have it back in service quickly without having to worry about where to store the remainder of the anesthetic you are using.

HIGHLY COMPATIBLE

Dräger vaporizers were designed to work with practically any Dräger anesthesia workstation. They come equipped with compatible connectors for both Dräger (DW-2000) and Selectatec (S-2000) as well as the standard 23mm ISO connector.

MADE WITH PATIENT SAFETY IN MIND

The color of each vaporizer corresponds to the color code of the anesthetic agent for which it was designed. The shape of each filling nozzle prevents filling the vaporizer with the wrong anesthetic agent. Our new Dräger AutoExclusion comfort

mounting also effectively eliminates the possibility of activating more than one vaporizer at a time. Once one vaporizer is turned on, the other is automatically blocked. Taken together, these features result in enhanced patient safety and a reduced likelihood of significant user error.

BACKUP POWER

The D-Vapor differs in function slightly. It is electrically driven and comes equipped with a battery which maintains function in case of a power failure.

SAFE AND SIMPLE TRANSPORT

When set in transport position, both the Dräger Vapor 2000 and the D-Vapor can withstand tipping or even being turned upside down without any ill effects. The anesthetic agent remains safely contained within the vaporizer. After connecting the vaporizer and releasing the transport position, they are immediately ready for operation according to specification.



D-877-2010

Dräger Vapor 2000 - isoflurane



D-876-2010


Dräger Vapor 2000 - sevoflurane




D-880-2010

Dräger D-Vapor - desflurane

D-90-2011

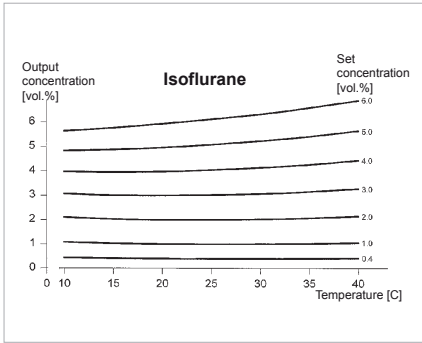

D-91-2011

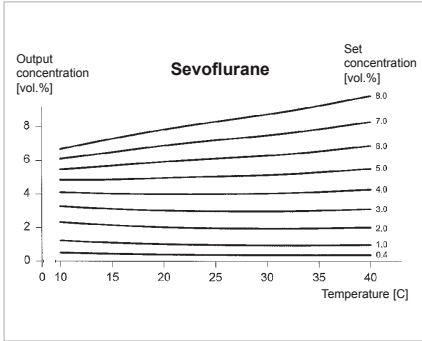

OPERATIONAL WITHIN A WIDER TEMPERATURE RANGE

The Dräger Vapor 2000 and D-Vapor can be safely operated at temperatures ranging from 10 – 40 °C. This makes them particularly suitable for use in operating rooms where cooling or heating is an essential factor, such as in pediatric surgery, cardiac surgery, burn units, and orthopedics.

ZERO CALIBRATION

The Vapor line of vaporizers requires absolutely no calibration: This saves you time.



D-7676-2010

APOLLO ANESTHESIA WORKSTATION RECORDING ANESTHETIC AGENTS

Apollo can automatically record total usage, patient uptake, and change of anesthetic agents during a case to be viewed at its conclusion.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
 Moislinger Allee 53–55
 23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Drägerwerk AG & Co. KGaA
 Moislinger Allee 53–55
 23558 Lübeck, Germany

USA

Dräger, Inc.
 3135 Quarry Road
 Telford, PA 18969-1042, USA
 Tel +1 800 4DRAGER
 (+1 800 437 2437)
 Fax +1 215 723 5935
info.usa@draeger.com

CANADA

Draeger Medical Canada Inc.
 2425 Skymark Avenue, Unit 1
 Mississauga, Ontario, L4W 4Y6
 Tel +1 905 212 6600
 Toll-free +1 866 343 2273
 Fax +1 905 212 6601
Canada.support@draeger.com

Locate your Regional Sales Representative at:
www.draeger.com/contact



Not all products, features, or services are for sale in all countries. Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.