

# Infinity® Explorer Software

Infinity® Explorer is a software application that lets you extend the viewing capability of Infinity monitors\* and integrate additional patient information on a separate display. Infinity Explorer runs on the Infinity C700 for IT workstation – the state-of-the-art medical grade computer component of the Infinity Omega solution.



Infinity Explorer software integrates real-time data from Infinity monitors and a wide variety of clinical applications at the point of care. The result? Faster and easier access to information, which supports better-informed clinical decision-making.

## BRINGS INFORMATION TO THE POINT OF CARE

Infinity Explorer software bridges the gap between patient monitors and the IT infrastructure of your hospital. The easy-to-navigate interface has up to 10 tabs – providing fingertip access to clinical applications. It lets you display all the parameters and waveforms that are available on the bedside monitor – including information from HL7 interfaces and data coming into the monitor from other point-of-care devices.

## SUPPORTS WORKFLOW AT THE ACUTE POINT OF CARE


Infinity Explorer enables dual control of monitoring functions from both the monitor and the workstation. You can control

alarms, display remote bedside monitors, and access the hospital's network – including PACS and clinical information systems – without the need for extra hardware.

## ENHANCED PATIENT MONITORING VIEWER

The Vitals screen allows interaction between a workstation and the connected bedside monitor, which is the primary monitoring device; offering the following information management capabilities:


- Event storage: 150 events, 20-second strips, all monitored waveforms
- Data storage: 96 hours of trended parameter information
- Data resolution: 30-second sampling
- Trend tables: 1, 5, 10, 15, 30 or 60-minute display formats
- Trend graphs: 1, 2, 8, 12, 24, 48, 72 or 96-hour display formats

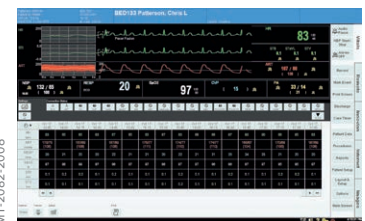

MT-2082-2008

### Vitals screen

Displays all monitored patient data and integrates it with information from other medical devices

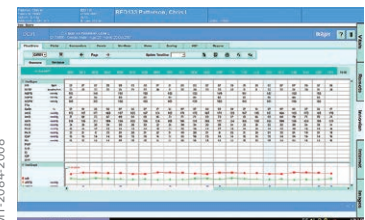

MT-2084-2008

### IT applications

Lets you view Web-based or Windows®-based applications of your choice



MT-2083-2008

### IT applications

Displays DICOM images and X-rays at the bedside

\* Infinity Omega includes a Delta or Delta XL monitor; Infinity Omega-S includes a Kappa monitor.

**SCALES TO MEET YOUR NEEDS**

Infinity Explorer software is designed as a set of applications, so you can begin with one application and then add others according to your needs.

**WORLD WIDE WEB VIEWER**

Lets you access the hospital network via Internet Explorer from the Vitals tab.

**REMOTE BEDSIDE VIEWER  
(REQUIRES INFINITY GATEWAY)**

Displays monitoring information from any remote bedside monitor that can be accessed by the Infinity Gateway server.

**WINDOWS/WEB-BASED CLINICAL APPLICATIONS**

Hosts up to 4 compatible Web-based applications, plus a maximum of 8 Windows-based applications.

**TERMINAL EMULATOR APPLICATIONS**

Acts as a terminal emulator for compatible hyper terminal applications.

**ORDERING INFORMATION**

Infinity® Explorer	MS18785
Omega	MS25200
Omega-S	MS25490
<b>Additional Software Applications</b>	
Remote Bedside Viewer (Requires Infinity® Gateway Server)	7488211
Windows/Web-based Clinical Applications	7488237
World Wide Web Viewer	7489292
Terminal Emulator Applications (VT52; VT100; ANSI; VTNT)	7488229
Video Option	MS17313

Infinity® is a registered trademark of Dräger.

Windows® is a registered trademark of Microsoft Corporation.

**CORPORATE HEADQUARTERS**

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany

[www.draeger.com](http://www.draeger.com)

**Manufacturer:**

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany

**REGION EUROPE CENTRAL  
AND EUROPE NORTH**

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany  
Tel +49 451 882 0  
Fax +49 451 882 2080  
[info@draeger.com](mailto:info@draeger.com)

**REGION EUROPE SOUTH**

Dräger Médical S.A.S.  
Parc de Haute Technologie  
d'Antony 2  
25, rue Georges Besse  
92182 Antony Cedex, France  
Tel +33 1 46 11 56 00  
Fax +33 1 40 96 97 20  
[dlmfr-contact@draeger.com](mailto:dlmfr-contact@draeger.com)

**REGION MIDDLE EAST, AFRICA**

Drägerwerk AG & Co. KGaA  
Branch Office  
P.O. Box 505108  
Dubai, United Arab Emirates  
Tel +971 4 4294 600  
Fax +971 4 4294 699  
[contactuae@draeger.com](mailto:contactuae@draeger.com)

**REGION ASIA / PACIFIC**

Draeger Medical  
South East Asia Pte Ltd.  
25 International Business Park  
#04-27/29 German Centre  
Singapore 609916, Singapore  
Tel +65 6572 4388  
Fax +65 6572 4399  
[asia.pacific@draeger.com](mailto:asia.pacific@draeger.com)

**REGION NORTH AMERICA**

Draeger Medical, Inc.  
3135 Quarry Road  
Telford, PA 18969-1042, USA  
Tel +1 215 721 5400  
Toll-free +1 800 437 2437  
Fax +1 215 723 5935  
[info.usa@draeger.com](mailto:info.usa@draeger.com)

**REGION CENTRAL  
AND SOUTH AMERICA**

Dräger Panama Comercial  
S. de R.L.  
Complejo Business Park,  
V tower, 10th floor  
Panama City  
Tel +507 377 9100  
Fax +507 377 9130  
[contactcsa@draeger.com](mailto:contactcsa@draeger.com)

Locate your Regional Sales  
Representative at:  
[www.draeger.com/contact](http://www.draeger.com/contact)

