


FitStar® NIV oronasal mask Consumables and Accessories

The FitStar disposable NIV mask has an ergonomic design and contains gel cushion technology that reduces pressure while maintaining an effective seal. Three color-coded elbow designs make choosing the right mask for CPAP/Bi-level, single or dual limb circuit NIV system easy.

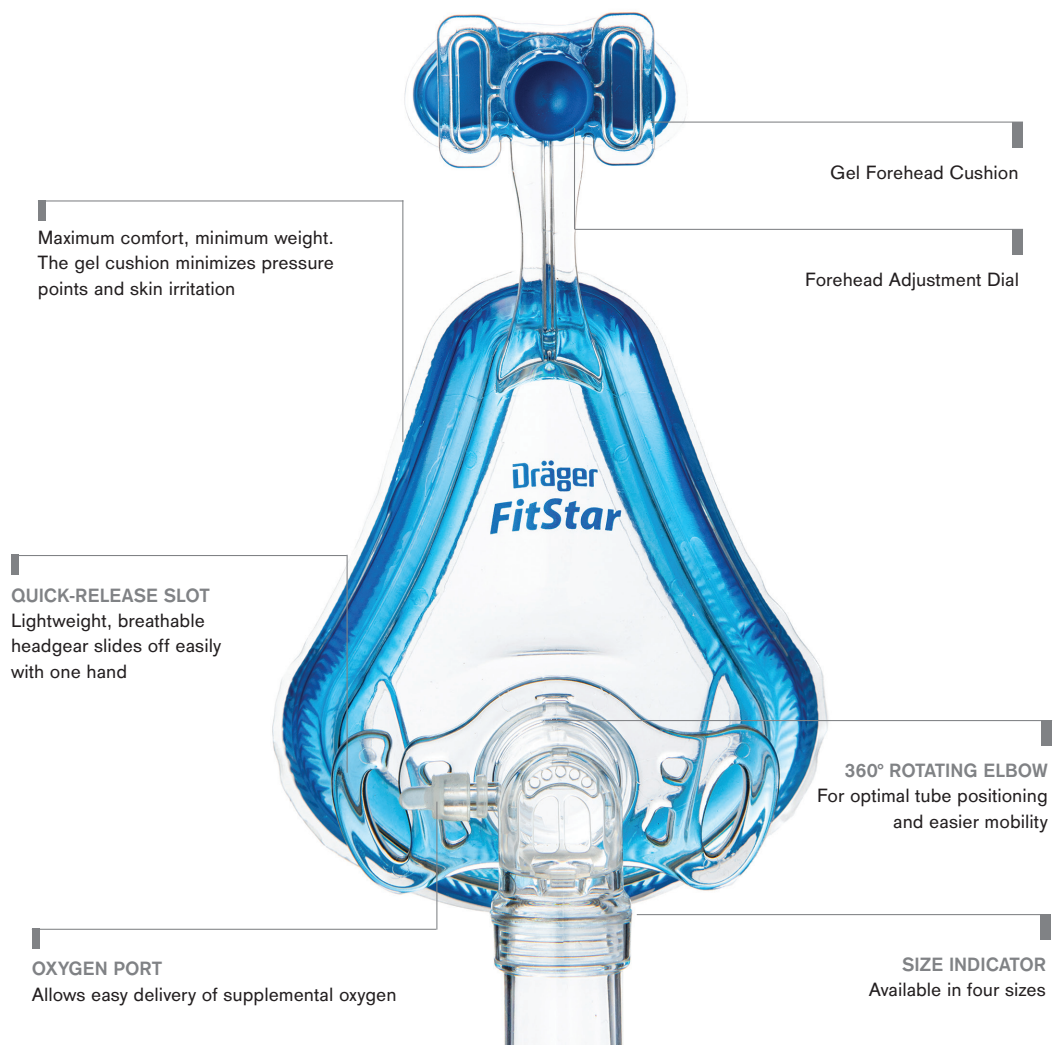

Benefits

Design

Designed with minimal intentional leaks and elbows that work with CPAP/Bi-level and single or dual limb circuits.

Does not contain natural rubber latex.

Adjustment Dial

Provides a wide, continuous range of adjustment ensuring a custom fit while removing pressure from the bridge of the nose.

Headgear

Maintains a secure and comfortable fit during treatment.

Size Indicator

Available in four sizes with simple and breathable quick-release headgear.

Details

A full line of disposable masks to meet your NIV needs, color-coded for easy identification.



Vented mask with an anti-asphyxia valve for use with CPAP or bi-level devices



Nonvented mask with an anti-asphyxia valve for use with single limb circuits



Nonvented mask for use with dual limb circuits

Ordering Information

FitStar masks with headgear specifications:

Product Number	Style	Leak Port	Entrianment Valve	Size	Mask Weight	Mask Dimensions
NA10700	Vented	✓	✓	Small	109.2g	3.60"D x 4.10"W x 7.25"H
NA10701	Vented	✓	✓	Medium	114.3g	3.60"D x 4.30"W x 7.40"H
NA10702	Vented	✓	✓	Large	119.4g	3.60"D x 4.38"W x 7.60"H
NA10703	Vented	✓	✓	Extra Large	120.1g	3.60"D x 4.42"W x 7.80"H
NA10704	Non-Vented	X	X	Small	108.6g	3.70"D x 4.10"W x 7.25"H
NA10705	Non-Vented	X	X	Medium	113.7g	3.70"D x 4.30"W x 7.40"H
NA10706	Non-Vented	X	X	Large	118.8g	3.70"D x 4.38"W x 7.60"H
NA10707	Non-Vented	X	X	Extra Large	119.5g	3.70"D x 4.42"W x 7.80"H
NA10708	Non-Vented	X	✓	Small	109.1g	3.65"D x 4.10"W x 7.25"H
NA10709	Non-Vented	X	✓	Medium	114.2g	3.65"D x 4.30"W x 7.40"H
NA10710	Non-Vented	X	✓	Large	119.3g	3.65"D x 4.38"W x 7.60"H
NA10711	Non-Vented	X	✓	Extra Large	120g	3.65"D x 4.42"W x 7.80"H

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Drägerwerk AG & Co. KGaA
Moislinger Allee 53-55
23542 Lübeck, Germany

USA

Draeger, Inc.
3135 Quarry Road
Telford, PA 18969-1042
Tel +1 800 4DRAGER
(+1 800 437 2437)
Fax +1 215 723 5935
info.usa@draeger.com