


O₂Star™ - Oxygen Therapy Accessories Consumables and Accessories

A comprehensive range of O₂ accessories providing quality and comfort for patients and caregivers alike.

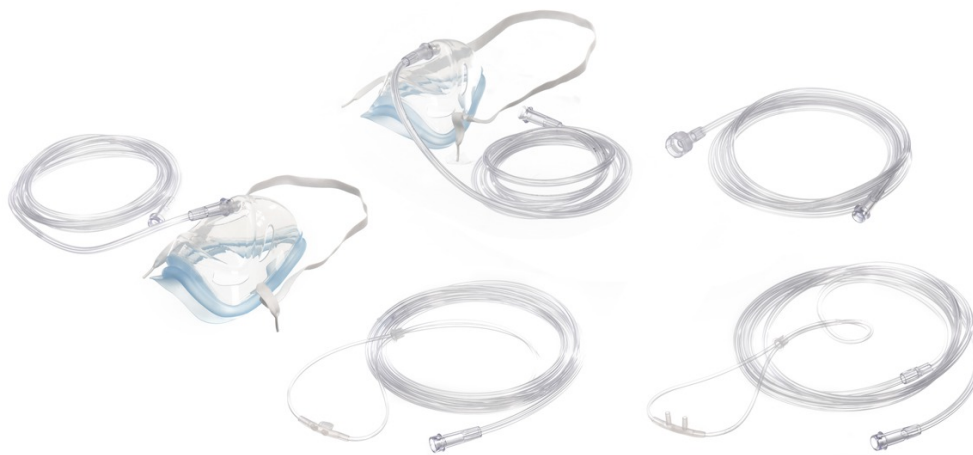

Benefits

Atmospheric oxygen

Room air contains approximately 21 % oxygen. This is generally enough for an individual to maintain adequate blood oxygenation through normal respiration.

Treatment of hypoxia

For various reasons, however, some medical conditions result in a shortage of oxygen to the tissues of the body (hypoxia). In many cases, this can be overcome by supplying oxygen to the patient at a higher concentration than it is found in the ambient air. In addition, the application of oxygen is often used during and following surgical procedures.

Specialized therapy

Oxygen can be delivered to the patient through tubing attached to a nasal cannula (prong), face mask, a face tent or a tracheostomy mask. Dräger provides a comprehensive range of oxygen delivery products for different applications, each particularly suited to treat a range of disorders resulting in hypoxia.

Nasal cannulas

Dräger's O₂Star nasal cannula range includes four different prong types including straight, curved, flared, and curved-flared shapes. Each is available in four different sizes ranging from small (S), medium (M), large (L) to extra-large (XL). The over-the-ear design maintains proper positioning of the soft nasal prong tips while allowing for complete freedom of movement. The attached oxygen tube is manufactured with a star-shaped lumen design; oxygen flow is maintained even if the tube is kinked.

Face masks

The shape of the O₂Star face mask allows for greater comfort and provides an effective seal. The light-weight design and in general wider, elastic head strap enable stable facial positioning. The head straps are reinforced at the ends, thus minimizing the risk of it being pulled through the mask. These face masks are available in two sizes - small (S) and medium/large (M/L).

Additional oxygen therapy products

Dräger's O₂Star product portfolio also includes Venturi mask kits with either six fixed Venturi valves or a single adjustable valve. Both Venturi mask kits are available in S and M/L sizes. Face tents and tracheostomy masks in two sizes (S and M/L) complete the product portfolio. Additional accessories such as oxygen tubes or various connectors are also available. They provide a secure connection to oxygen tubes or breathing circuits.

Benefits

Color-coded, light-weight

The O₂Star product portfolio continues with the “traffic-light” color coding principle – well-known from other disposable Dräger mask lines. The product labels are color-coded; smaller sizes have a red label, medium sizes a yellow label, and larger sizes a green label. This approach allows for easy size identification, even from a distance. The light-weight, transparent packaging bag is marked and designed for fast and easy opening. All materials are heat resistant and the complete portfolio is entirely phthalate-free. The oxygen face masks are completely PVC-free. None of the O₂Star products contains natural rubber.

Ordering Information

| DESCRIPTION | PACKAGING UNIT | ORDER NO. |
|-----------------------------------------------------------------------------------------------|----------------|-----------|
| O ₂ Star™ Non-rebreather oxygen mask with connection tube 2.1 m, disposable, M/L | 50 | MP01920 |
| O ₂ Star™ Medium concentration oxygen mask with conn. tube 2.1 m, disposable, M/L | 50 | MP01921 |
| O ₂ Star™ Venturi mask kit, adjustable venturi valve, disposable, M/L | 25 | MP01922 |
| O ₂ Star™ Venturi mask kit, fixed venturi valve, disposable, M/L | 25 | MP01923 |
| O ₂ Star™ Tracheostomy mask, disposable, M/L | 50 | MP01924 |
| O ₂ Star™ Face tent, disposable, M/L | 50 | MP01925 |
| O ₂ Star™ Oxygen and aerosol mask with 22 mm OD, disposable, M/L | 50 | MP01926 |
| O ₂ Star™ Oxygen connection tube, 2.1 m, disposable | 25 | MP01927 |
| O ₂ Star™ Venturi valve kit, 6 fixed venturi valves, disposable | 5 | MP01928 |
| O ₂ Star™ Adjustable venturi valve, disposable | 5 | MP01929 |
| O ₂ Star™ Non-rebreather oxygen mask with connection tube 2.1 m, disposable, S | 50 | MP01930 |
| O ₂ Star™ Medium concentration oxygen mask w connection tube 2.1 m, disposable, S | 50 | MP01931 |
| O ₂ Star™ Tracheostomy mask, disposable, S | 50 | MP01932 |
| O ₂ Star™ Face tent, disposable, S | 50 | MP01933 |
| O ₂ Star™ Oxygen and aerosol mask with 22 mm OD, disposable, S | 50 | MP01934 |
| O ₂ Star™ Venturi mask kit, adjustable venturi valve, disposable, S | 25 | MP01935 |
| O ₂ Star™ Venturi mask kit, fixed venturi valve, disposable, S | 25 | MP01936 |
| O ₂ Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, S | 50 | MP01920 |
| O ₂ Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, M | 50 | MP02921 |
| O ₂ Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, L | 50 | MP02922 |
| O ₂ Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, XL | 50 | MP02923 |
| O ₂ Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, S | 50 | MP02924 |
| O ₂ Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, M | 50 | MP02925 |
| O ₂ Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, L | 50 | MP02926 |

| | | |
|------------------------------------------------------------------------------------------------------|----|---------|
| O ₂ Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, XL | 50 | MP02927 |
| O ₂ Star™ Oxygen nasal cannula flared with connection tube 2.1 m, disposable, S | 50 | MP02928 |
| O ₂ Star™ Oxygen nasal cannula flared with connection tube 2.1 m, disposable, M | 50 | MP02929 |
| O ₂ Star™ Oxygen nasal cannula flared with connection tube 2.1 m, disposable, L | 50 | MP02930 |
| O ₂ Star™ Oxygen nasal cannula flared with connection tube 2.1 m, disposable, XL | 50 | MP02931 |
| O ₂ Star™ Oxygen nasal cannula curved-flared with conn. tube 2.1 m, disposable, S | 50 | MP02932 |
| O ₂ Star™ Oxygen nasal cannula curved-flared with conn. tube 2.1 m, disposable, M | 50 | MP02933 |
| O ₂ Star™ Oxygen nasal cannula curved-flared with conn. tube 2.1 m, disposable, L | 50 | MP02934 |
| O ₂ Star™ Oxygen nasal cannula curved-flared with conn. tube 2.1 m, disposable, XL | 50 | MP02935 |
| O ₂ Star™ Connector kit: 6mm OD/22mm OD + 6mm OD/22mm ID with connection tube, disposable | 5 | MP02940 |
| O ₂ Star™ Connector kit for breathing circuit for Aquapor, disposable | 5 | MP02941 |
| O ₂ Star™ Connector 22 mm ID / 22 mm ID, disposable | 5 | MP02942 |
| O ₂ Star™ Oxygen connection tube, 4 m, disposable | 25 | MP01943 |
| O ₂ Star™ Oxygen connection tube 2.1 m with 9/16" connector, disposable | 25 | MP01944 |
| O ₂ Star™ Oxygen connection tube 4 m with 9/16" connector, disposable | 25 | MP01945 |
| O ₂ Star™ Demo package oxygen therapy accessories incl. mask, nasal cannula, con. tube | 5 | MP01946 |

Notes

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
www.draeger.com

Manufacturer:

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23542 Lübeck, Germany

REGION DACH

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA PACIFIC

Draeger Singapore Pte. Ltd.
25 International Business Park
#04-20/21 German Centre
Singapore 609916
Tel +65 6308 9400
Fax +65 6308 9401
asia.pacific@draeger.com

REGION CENTRAL AND SOUTH AMERICA

Dräger Panama S. de R.L.
59 East Street, Nuevo Paitilla
House 30, San Francisco Town
Panama City, Panama
Tel +507 377 9100
Fax +507 377 9130
servicioalcliente@draeger.com

Locate your Regional Sales
Representative at:
www.draeger.com/contact

